

3dDeconvolve

Advanced Features

**Just in case you weren't
confused enough already**

Other Features of 3dDeconvolve - 1

- **-input1D** = used to process a single time series, rather than a dataset full of time series
 - ★ e.g., test out a stimulus timing sequence on sample data
 - ★ **-nodata** option can be used to check for collinearity
- **-censor** = used to turn off processing for some time points
 - ★ for time points that are “bad” (e.g., too much movement; scanner hiccup)
 - ★ **-CENSORTR 2:37** = newer way to specify omissions (e.g., run #2, index #37)
- **-sresp** = output standard deviation of HRF (β) estimates
 - ★ can then plot error bands around HRF in AFNI graph viewer
- **-errts** = output residuals (difference between fitted model and data)
 - ★ for statistical analysis of time series noise
- **-TR_times dt** = calculate **-iresp** and **-sresp** HRF results with time step **dt** (instead of input dataset TR)
 - ★ Can be used to make HRF graphs look better
- **-jobs N** = run with independent threads — **N** of them
 - ★ extra speed, if you have a dual-CPU system (or more)!

Other Features - 2

<http://afni.nimh.nih.gov/pub/dist/doc/misc/Decon/DeconSummer2004.html>

<http://afni.nimh.nih.gov/pub/dist/doc/misc/Decon/DeconSpring2007.html>

- Equation solver: Program computes **condition number** for **X** matrix (measures of how sensitive regression results are to changes in **X**)
 - ★ If the condition number is “bad” (too big), then the program will not actually proceed to compute the results
 - ★ You can use the **-GOFORIT** option on the command line to force the program to run despite **X** matrix warnings
 - But you should strive to understand why you are getting these warnings!!
- Other matrix checks:
 - ★ Duplicate stimulus filenames, duplicate regression matrix columns, all zero matrix columns
- ★ Check the screen output for **WARNINGS**s and **ERRORS**★
 - ★ Such messages also saved into file **3dDeconvolve.err**

Other Features - 3

- All-zero regressors *are* allowed (via `-allzero_OK` or `-GOFORIT`)
 - ★ Will get zero weight in the solution
 - ★ Example: task where subject makes a choice for each stimulus (e.g., male or female face?)
 - You want to analyze correct and incorrect trials as separate cases
 - What if some subject makes no mistakes? Hmmm...
 - Can keep the all-zero regressor (e.g., all `-stim_times = *`)
 - Input files and output datasets for error-making and perfect-performing subjects will be organized the same way
-
- **3dDeconvolve_f** program can be used to compute linear regression results in single precision (7 decimal places) rather than double precision (16 places)
 - ★ For better speed, but with lower numerical accuracy
 - ★ Best to do at least one run ***both*** ways to check if results differ significantly (Equation solver *should* be safe, but ...)

Other Features - 4


- Default output format is 16-bit short integers, with a scaling factor for each sub-brick to convert it to floating point values
 - ★ `-float` option can be used to get 32-bit floating point format output — more precision, and more disk space

- `3dDeconvolve` recommends a `-polort` value, and prints that out as well as the value you chose (or defaulted to)
 - ★ `-polort A` can be used to let the program set the detrending (AKA “high pass filtering”, since detrending removes low frequency content from data) level automatically

- `-stim_file` is used to input a column directly into `X` matrix
 - ★ Motion parameters (as in previous examples)
 - ★ If you create a stimulus+response model outside `3dDeconvolve` (e.g., using program `waver`)

Other Features - 5

- **-stim_times** has some other basis function options for the HRF model besides **BLOCK** and **TENT**
 - ★ **CSPLIN** = cubic spline instead of **TENT** = linear spline
 - Same parameters: (**start, stop, number of regressors**)
 - Can be used as a “drop in” replacement for **TENT**

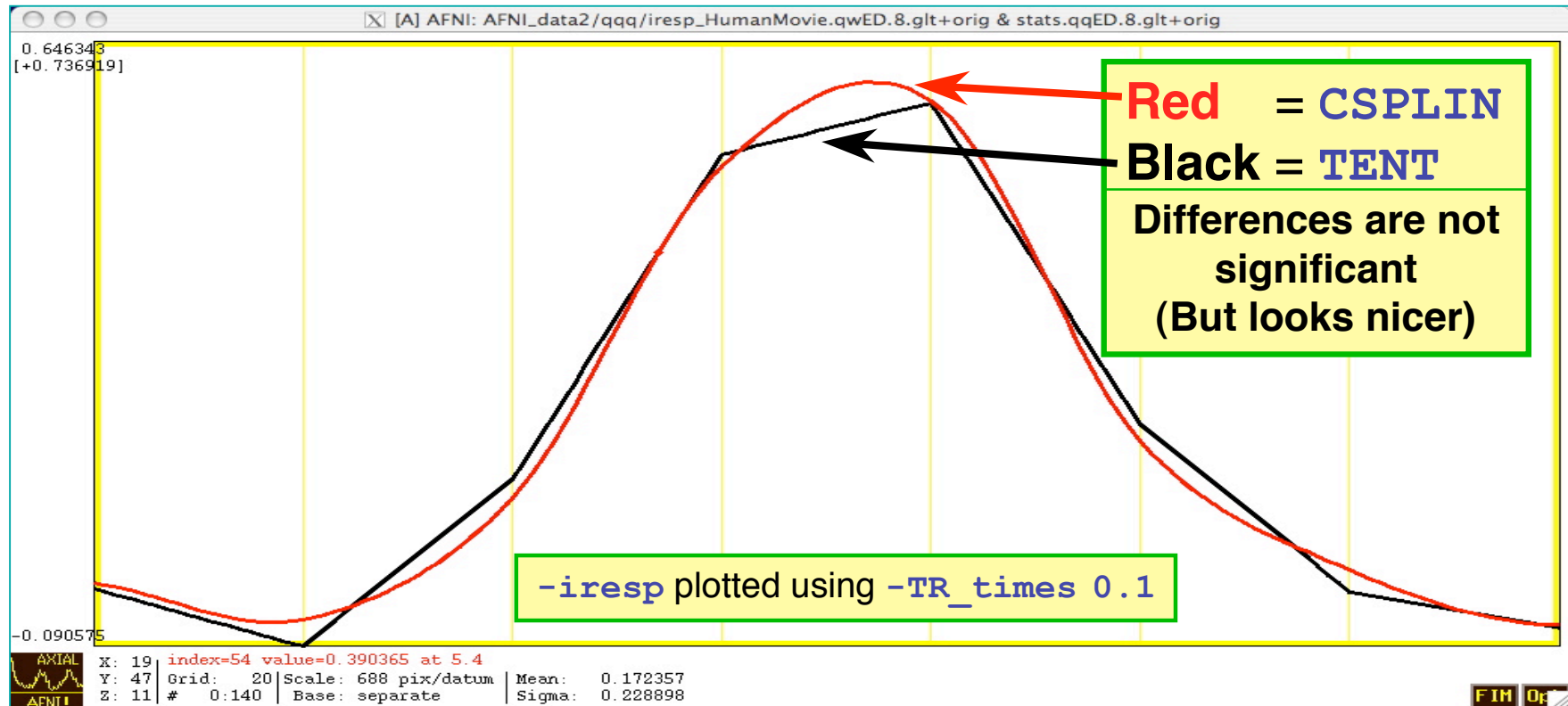

Other Features - 6

- **-fitts** option is used to create a synthetic dataset
 - ★ each voxel time series is full (signal+baseline) model as fitted to the data time series in the corresponding voxel location
- **3dSynthesize** program can be used to create synthetic datasets from *subsets* of the full model
 - ★ Uses **-x1D** and **-cbucket** outputs from **3dDeconvolve**
 - **-cbucket** stores β coefficients for each **X** matrix column into dataset
 - **-x1D** stores the matrix columns (and **-stim_labels**)
 - ★ Potential uses:
 - Baseline only dataset
 - ➔ **3dSynthesize -cbucket fred+orig -matrix fred.x1D -select baseline -prefix fred_base**
 - ➔ Could subtract this dataset from original data (via **3dcalc**) to get signal+noise dataset that has no baseline component left
 - Just one stimulus class model (+ baseline) dataset
 - ➔ **3dSynthesize -cbucket fred+orig -matrix fred.x1D -select baseline Faces -prefix fred_Faces**

Other Recent Small Changes

- Defaults are changed:
 - ★ **-nobout** & **-full_first** & **-bucket** & **-x1D** are always implied
 - ★ Names of statistics sub-bricks are slightly altered (to be more consistent)

- Checks if **-stim_times** inputs are out of range (AKA: the PSFB syndrome)
 - ★ Prints **WARNING** message, but continues analysis

- When using **-nodata** with **-stim_times**, it is important to give the number of time points and the TR, as in **-nodata 250 2.3**
 - ★ With **-input1D**, use **-TR_1D 2.3** to specify TR


IM Regression - 1

- **IM** = Individual **M**odulation
 - ★ Compute *separate* amplitude of response for each stimulus
 - Instead of computing average amplitude of responses to multiple stimuli in the same class
 - ★ Response amplitudes (β_s) for each individual block/event will be highly noisy
 - Can't use individual activation map for much
 - Must pool the computed β_s in some further statistical analysis (t -test via **3dttest**? inter-voxel correlations in the β_s ? correlate β_s with something else?)
 - ★ Usage: **-stim_times_IM k tname model**
 - Like **-stim_times**, but creates a separate regression matrix column for each time given

IM Regression - 2

- Only application of IM thus far has been in checking some data we received from another institution
- Experiment: 64 blocks of sensorimotor task (8 runs each with 8 blocks)

Plot of 64 **BLOCK** β s from **-cbucket** output

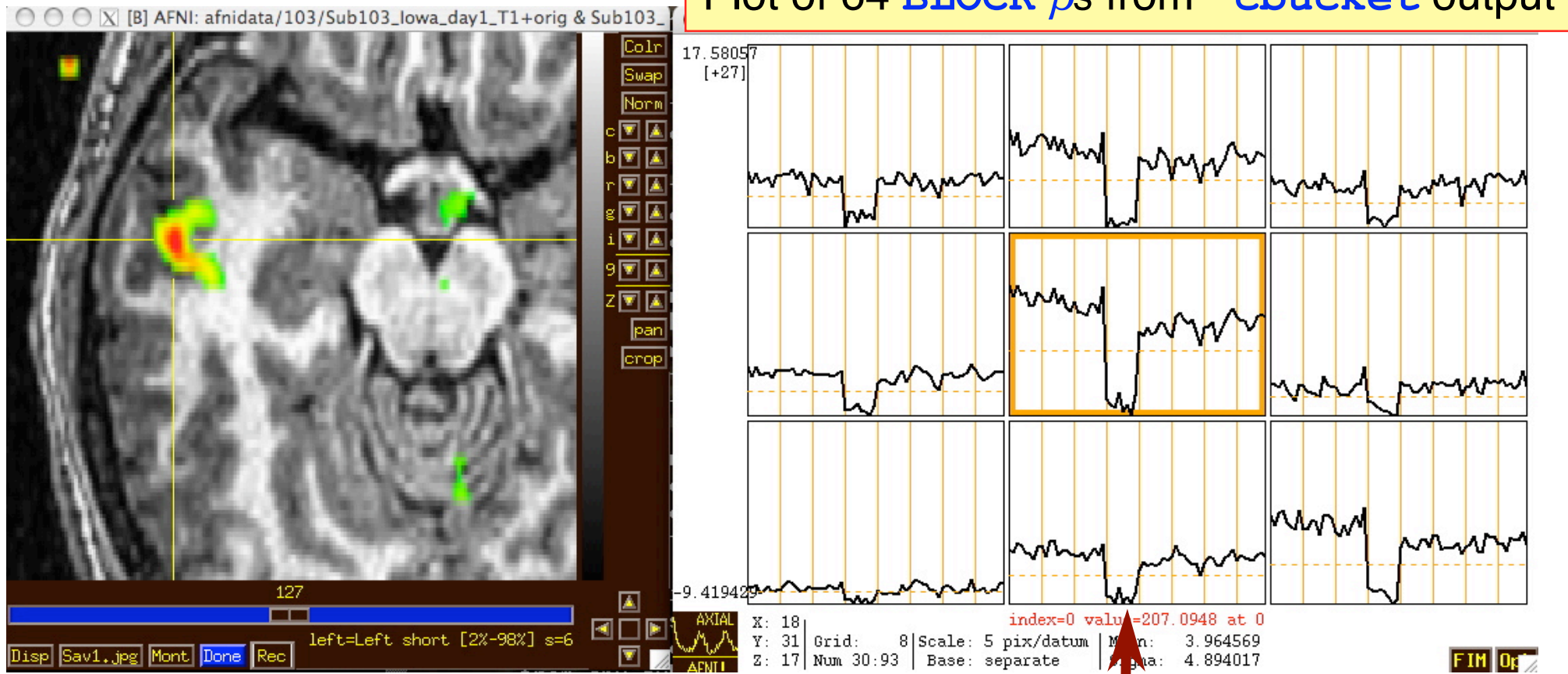

N.B.: sign reversal in run #4 = stimulus timing error!

IM Regression - 3

- IM works naturally with blocks, which only have 1 amplitude parameter per stimulus
- With event-related experiment and deconvolution, have multiple amplitude parameters per stimulus
 - ★ Difficulty: each event in same class won't get the same shaped HRF this way
 - ★ Desideratum: allow response shape to vary (that's deconvolution), but only allow amplitude to vary between responses in the same stimulus class
 - ★ Problem: get unknowns that multiply each other (shape parameters \times amplitude parameters) — and we step outside the realm of *linear* analysis
 - ★ Possible solution: **semi-linear** regression (nonlinear in global shape parameters, linear in local amplitude params)

AM Regression - 1

- **AM** = **A**mplitude **M**odulated (or **M**odulation)
 - ★ Have some extra data measured about each response to a stimulus, and *maybe* the BOLD response amplitude is modulated by this
 - ★ Reaction time; Galvanic skin response; Pain level perception; Emotional valence (happy or sad or angry face?)
 - Want to see if some brain activations vary proportionally to this **ABI** (**A**uxiliary **B**ehavioral **I**nformation)
-
- Discrete levels (2 or maybe 3) of ABI:
 - ★ Separate the stimuli into sub-classes that are determined by the ABI (“on” and “off”, maybe?)
 - ★ Use a GLT to test if there is a difference between the fMRI responses in the sub-classes

```
3dDeconvolve ... \
  -stim_times 1 regressor_on.1D 'BLOCK(2,1)' -stim_label 1 'On' \
  -stim_times 2 regressor_off.1D 'BLOCK(2,1)' -stim_label 2 'Off' \
  -gltsym 'SYM: +On | +Off' -glt_label 1 'On+Off' \
  -gltsym 'SYM: +On -Off' -glt_label 2 'On-Off' ...
```

- “**On+Off**” tests for any activation in *either* the “on” or “off” conditions
- “**On-Off**” tests for differences in activation *between* “on” and “off” conditions
- Can use **3dcalc** to threshold on *both* statistics at once to find a **conjunction**

AM Regression - 2

- Continuous (or several finely graded) ABI levels
 - ★ Want to find active voxels whose activation level also depends on ABI
 - ★ **3dDeconvolve** is a linear program, so must make the assumption that the change in FMRI signal as ABI changes is linearly proportional to the changes in the ABI values
- Need to make 2 separate regressors
 - ★ One to find the mean FMRI response (the usual `-stim_times` analysis)
 - ★ One to find the variations in the FMRI response as the ABI data varies
- The second regressor should have the form

$$r_{AM2}(t) = \sum_{k=1}^K h(t - \tau_k) \cdot (a_k - \bar{a})$$

- ★ Where a_k = value of k^{th} ABI value, and \bar{a} is the average ABI value
- Response (β) for first regressor is standard activation map
- Statistics and β for second regressor make activation map of places whose BOLD response changes with changes in ABI
 - ★ Using 2 regressors allows separation of voxels that are active but are *not* detectably modulated by the ABI from voxels that *are* ABI-sensitive

AM Regression - 3

- New feature of **3dDeconvolve**: `-stim_times_AM2`
- Use is very similar to standard `-stim_times`
 - ★ `-stim_times_AM2 1 times_ABI.1D 'BLOCK(2,1)'`
 - ★ The `times_ABI.1D` file has time entries that are “married” to ABI values:


```
10*5  23*4  27*2  39*5
17*2  32*5
*
16*2  24*3  37*5  41*4
```
 - ★ Such files can be created from 2 standard ASCII .1D files using the new **1dMarry** program
 - The `-divorce` option can be used to split them up
- **3dDeconvolve** automatically creates the two regressors (unmodulated and amplitude modulated)
 - ★ Use `-fout` option to get statistics for activation of the pair of regressors (i.e., testing null hypothesis that *both* β weights are zero: that there is no ABI-independent *or* ABI-proportional signal change)
 - ★ Use `-tout` option to test each β weight separately
 - ★ Can **1dplot** **X** matrix columns to see each regressor

AM Regression - 4

- The **AM** feature is new, and so needs some practical user experiences before it can be considered “standard practice”
 - ★ In particular: don’t know how much data or how many events are needed to get good ABI-dependent statistics
- If you want, **-stim_times_AM1** is also available
 - ★ It only builds the regressor proportional to ABI data directly, with no mean removed:


$$r_{\text{AM1}}(t) = \sum_{k=1}^K h(t - \tau_k) \cdot a_k$$
 - ★ Can’t imagine what value this option has, but you never know ... (if you can think of a good use, let me know)
- Future directions:
 - ★ Allow more than one amplitude to be married to each stimulus time (insert obligatory polygamy/polyandry joke here)
 - How many ABI types at once is too many? I don’t know.
 - ★ How to deal with unknown nonlinearities in the BOLD response to ABI values? I don’t know. (Regress each event separately, then compute MI?)
 - ★ Deconvolution with amplitude modulation? Requires more thought.

AM Regression - 5

Timing: AM.1D = 10*1 30*2 50*3 70*1 90*2 110*3 130*2 150*1 170*2 190*3 210*2 230*1


- 3dDeconvolve -nodata 300 1.0 -num_stimts 1 \
 - stim_times_AM1 1 AM.1D 'BLOCK(10,1)' -x1D AM1.x1D
- 1dplot AM1.x1D'[2]'

AM1 model of signal
(modulation = ABI)

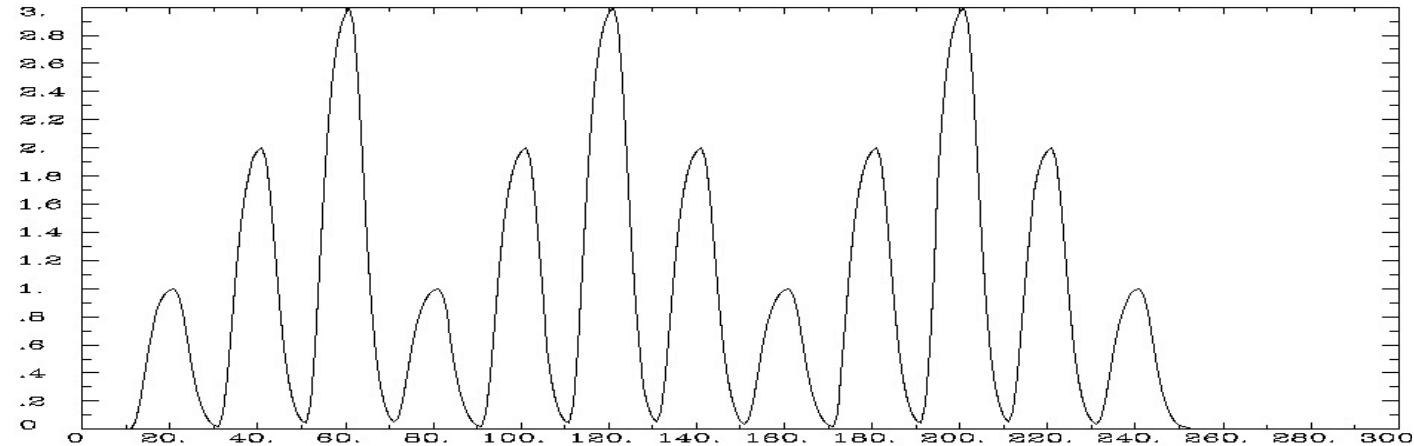



- 3dDeconvolve -nodata 300 1.0 \
 - num_stimts 1 \
 - stim_times_AM2 1 \
 - AM.1D 'BLOCK(10,1)' \
 - x1D AM2.x1D
- 1dplot -sepscl \
 - AM2.x1D'[2,3]'

AM2 model of signal:
is 2D sub-space
spanned by these 2
time series

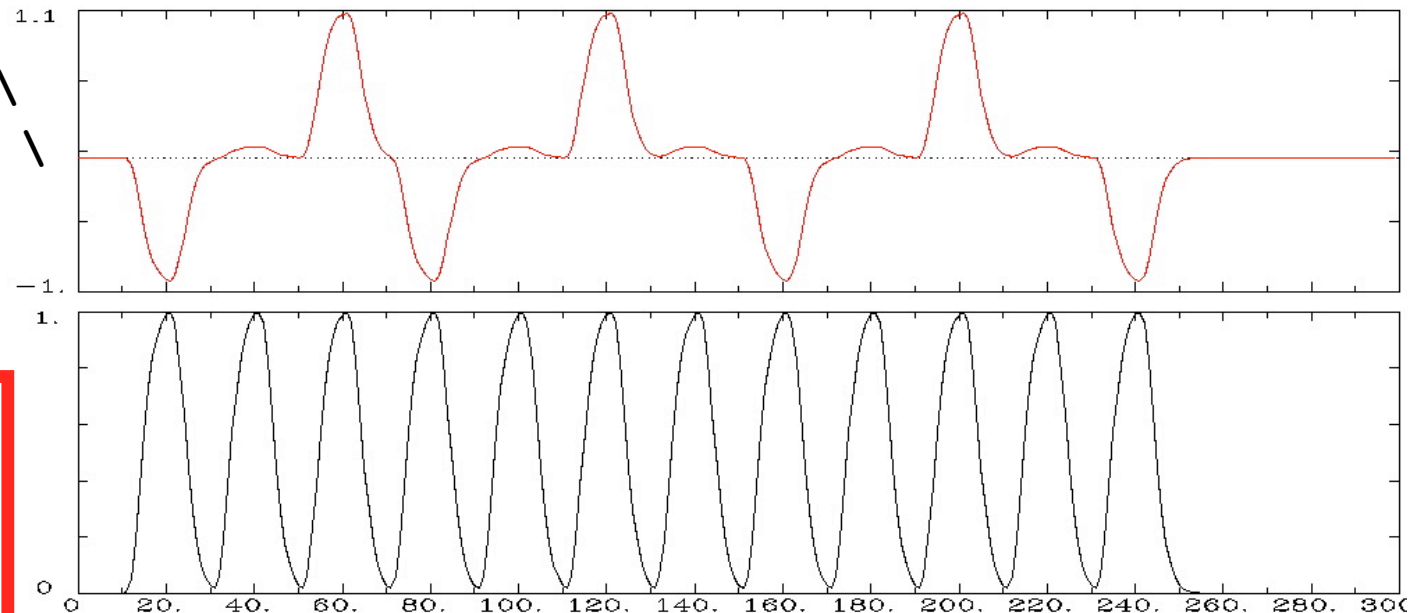

AM Regression - 6


- First actual user: Whitney Postman (formerly NIDCD; PI=Al Braun)
- Picture naming task in aphasic stroke patient
- ABI data = number of alternative names for each image (e.g., “balcony” & “porch” & “veranda”, vs. “strawberry”), from 1 to 18
 - 8 imaging runs, 144 stimulus events
- 2 slices showing activation map for BOLD responses proportional to ABI (β_{AM2})
 - What does this mean? Don't ask me!

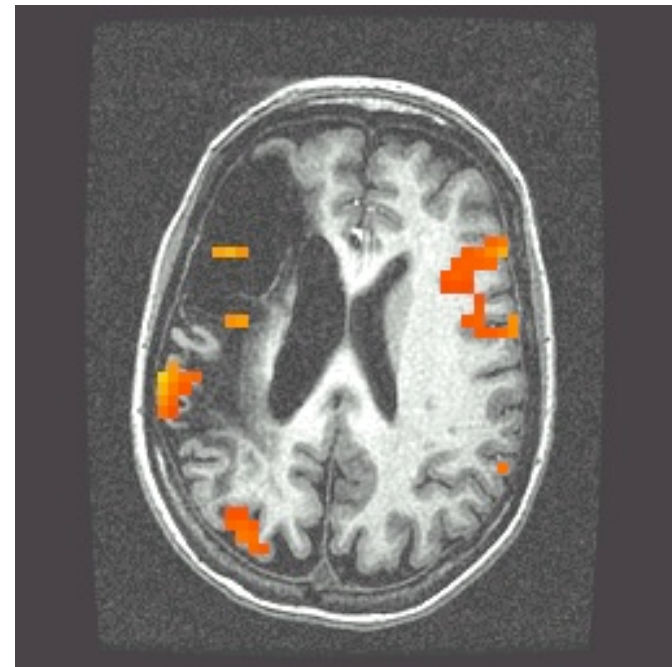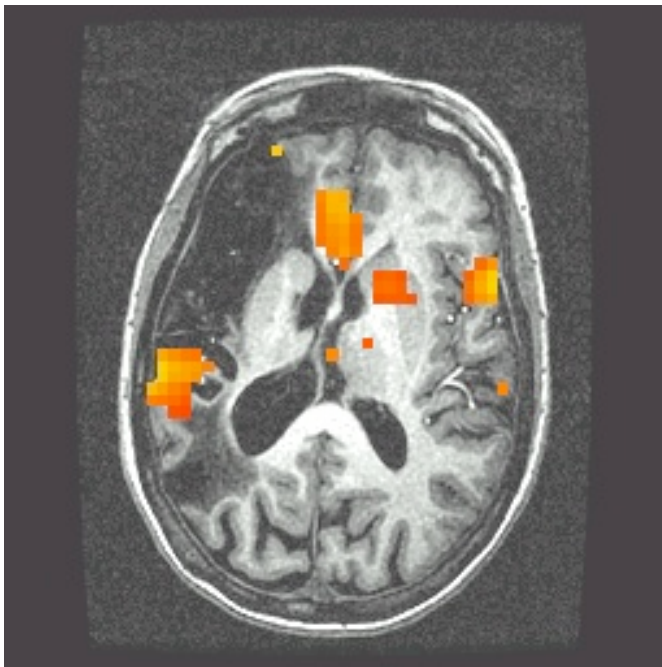


AM Regression - 7

- Alternative: use **IM** to get individual β s for each block/event and then do external regression statistics on those values
- Could do nonlinear fitting via **3dNIfim**, or inter-class contrasts via **3dttest**, **3dLME**, **3dANOVA**, etc.
- What is better: **AM** or **IM++** ?
 - We don't know – experience with these options is limited thus far – you can always try both!
 - If **AM** doesn't fit your models/ideas, then **IM** is clearly the way to go
 - Probably need to consult with SSCC to get some hints/advice

Other Advanced Topics in Regression

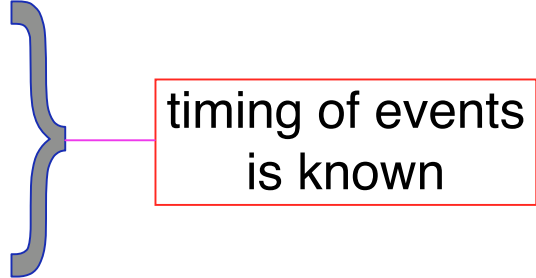- Can have activations with multiple phases that are not always in the same time relationship to each other; e.g.:
 - a) subject gets cue #1
 - b) variable waiting time (“hold”)
 - c) subject gets cue #2, emits response
 - which depends on both cue #1 and #2
- 
- timing of events is known
- ★ Cannot treat this as one event with one HRF, since the different waiting times will result in different overlaps in separate responses from cue #1 and cue #2
 - ★ Solution is multiple HRFs: separate HRF (fixed shape or deconvolution) for cue #1 times and for cue #2 times
 - Must have significant variability in inter-cue waiting times, or will get a nearly-collinear model
 - impossible to tell tail end of HRF #1 from the start of HRF #2, if always locked together in same temporal relationship
 - How much variability is “significant”? Good question.

Even More Complicated Case

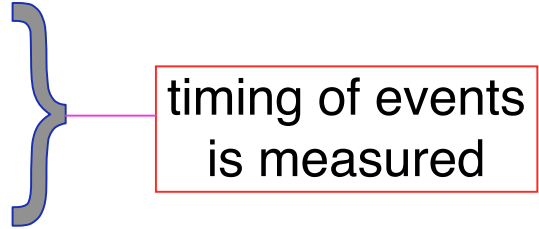- Solving a visually presented puzzle:
 - a) subject sees puzzle
 - b) subject cogitates a while
 - c) subject responds with solution

timing of events
is measured
- The problem is that we expect some voxels to be significant in phase (b) as well as phases (a) and/or (c)
- Variable length of phase (b) means that shape for its response varies between trials
 - ★ Which is contrary to the whole idea of averaging trials together to get decent statistics (which is basically what linear regression for the β weights does, in an elaborate sort of way)
- Could assume response **amplitude** in phase (b) is constant across trials, and response **duration** varies directly with time between phases (a) and (c)
 - ★ Need three HRFs
 - ★ Can't generate (b) HRF in **3dDeconvolve**

Noise Issues

- “Noise” in fMRI is caused by several factors, not completely characterized
 - ★ MR thermal noise (well understood, unremovable)
 - ★ Cardiac and respiratory cycles (partly understood)
 - In principle, could measure these sources of noise separately and then try to regress them out
 - ➔ RETROICOR program underway (Rasmus Birn of FIM/NIMH)
 - ★ Scanner fluctuations (e.g., thermal drift of hardware)
 - ★ Small subject head movements (10-100 mm)
 - ★ Very low frequency fluctuations (periods longer than 100 s)
- Data analysis should try to remove what can be removed and allow for the statistical effects of what can't be removed
 - ★ “Serial correlation” in the noise time series affects the t - and F -statistics calculated by **3dDeconvolve**
 - ★ At present, nothing is done to correct for this effect (by us)

Nonlinear Regression

- Linear models aren't the only possibility
 - ★ e.g., could try to fit HRF of the form $h(t) = a \cdot t^b \cdot e^{-t/c}$
 - ★ Unknowns b and c appear nonlinearly in this formula
- Program **3dNLFim** can do nonlinear regression (including nonlinear deconvolution)
 - ★ User must provide a C function that computes the model time series, given a set of parameters (e.g., a , b , c)
 - We could help you develop this C model function
 - Several sample model functions in the AFNI source code distribution
 - ★ Program then drives this C function repeatedly, searching for the set of parameters that best fit each voxel
 - ★ Has been used to fit pharmacological wash-in/wash-out models (difference of two exponentials) to fMRI data acquired during pharmacological challenges
 - e.g., injection of nicotine, cocaine, ethanol, etc.
 - these are difficult experiments to do **and** to analyze

Spatial Models of Activation

- Smooth data in space before analysis

- Average data across anatomically-selected regions of interest ROI (before or after analysis)
 - Labor intensive (*i.e.*, hire more students)

- Reject isolated small clusters of above-threshold voxels after analysis

Spatial Smoothing of Data

- Reduces number of comparisons
- Reduces noise (by averaging)
- Reduces spatial resolution
 - Blur it enough: Can make fMRI results look like low resolution PET data
- Smart smoothing: average **only** over nearby brain or gray matter voxels
 - Uses resolution of fMRI cleverly
 - New AFNI program: **3dBlurToFWHM**
 - Or: average over selected ROIs
 - Or: cortical surface based smoothing

Spatial Clustering


- Analyze data, create statistical map (e.g., t statistic in each voxel)

- Threshold map at a low t value, in each voxel separately
 - Will have many false positives

- Threshold map by rejecting clusters of voxels below a given size

- Can control false-positive rate by adjusting t threshold and cluster-size thresholds together

Cluster-Based Detection

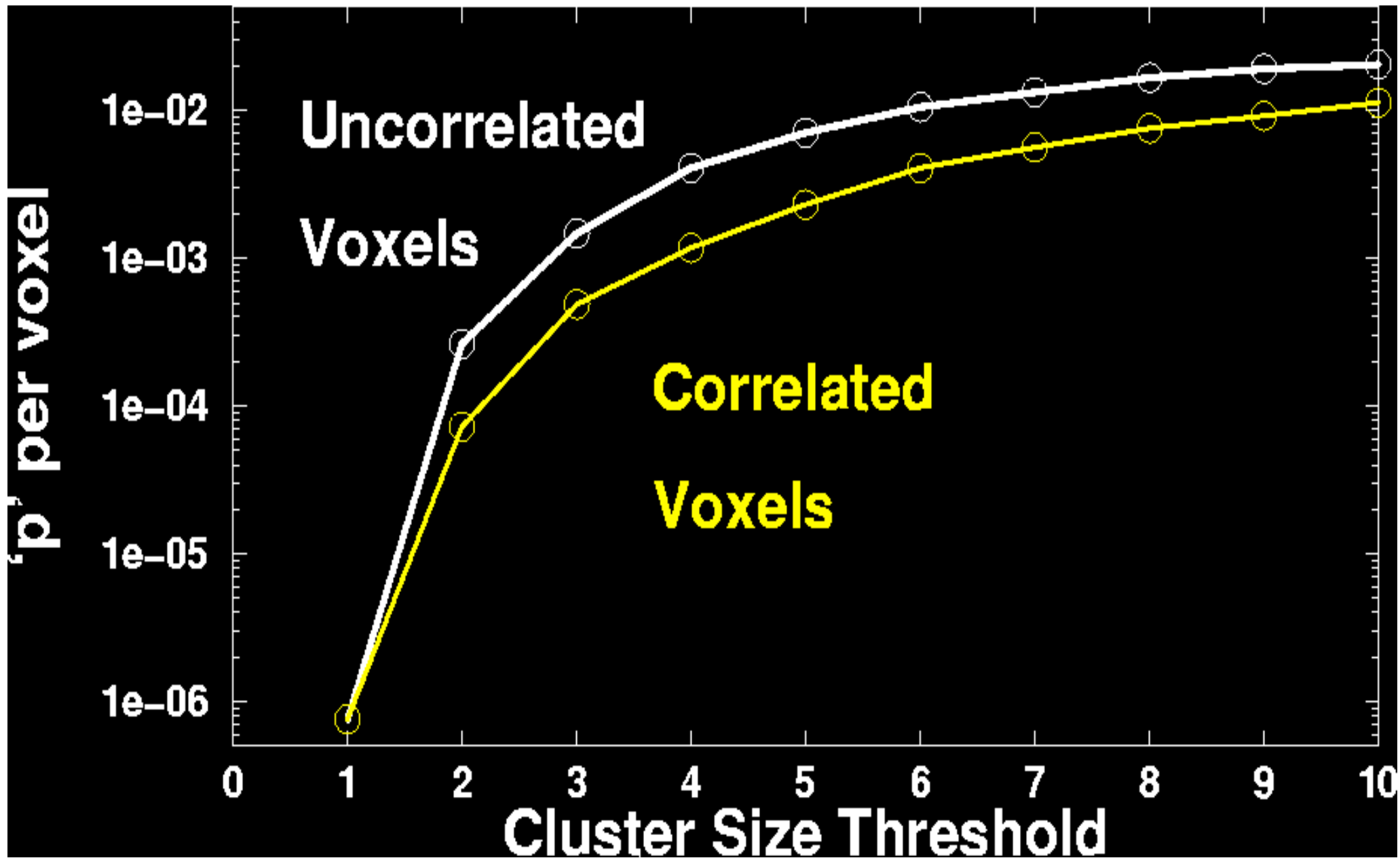


What the World Needs Now

- Unified HRF/Deconvolution ⊕ Blob analysis
 - Time ⊕ Space patterns computed all at once, instead of arbitrary spatial smoothing
 - Increase statistical power by bringing data from multiple voxels together cleverly
 - Instead of time analysis followed by spatial analysis (described earlier)
 - Instead of component-style analyses (e.g., ICA) that do not use stimulus timing
- Difficulty: models for spatial blobs
 - Little information *à priori* ⇒ must be adaptive

3dBlurToFWHM

- New program to smooth FMRI time series datasets to a specified smoothness (as estimated by FWHM of noise spatial correlation function)
 - ★ Don't just add smoothness (à la **3dmerge**) but control it (locally and globally)
 - ★ Goal: use datasets from diverse scanners
- Why blur FMRI time series?
 - ★ Averaging neighbors will reduce noise
 - ★ Activations are (usually) blob-ish (several voxels across)
 - ★ Diminishes the multiple comparisons problem
- **3dBlurToFWHM** blurs only inside a mask
 - ★ To avoid mixing air (noise-only) and brain voxels
 - ★ Partial Differential Equation (PDE) based blurring method
 - 2D (intra-slice) or 3D blurring

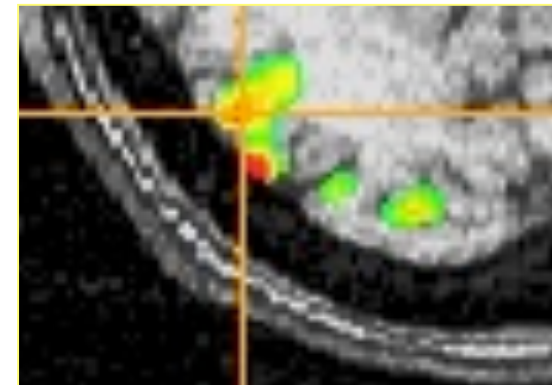

In the Pondering Stages

- “Area under curve” addition to `-gltsym` to allow testing of pieces of HRF models from `-stim_times`
- Slice- and/or voxel-dependent regressors
 - ★ For physiological noise cancellation, etc.
 - ★ To save memory? (Could process each slice separately)
 - One slice-at-a-time regression can be done in a Unix script, using 3dZcutup and 3dZcat programs
- Extend AM regression to allow for more than 1 piece of auxiliary information at each stimulus time
- Interactive tool to examine `-x1D` matrix for problems
 - ★ and `3dDeconvolve` testing of GLT submatrices
- Semi-linear deconvolution program