


Creating Datasets with Program to3d

- to3d reads image files — each containing 1 or more 2D slices — and assembles them into AFNI datasets
- The collection of all the 2D slice data forms the .BRIK file
 - an AFNI dataset must have at least 2 slices!
 - if you have single slice data, you could duplicate the slices
- You must also provide to3d with some auxiliary data (for the .HEAD file):
 - ◇ Orientation of the slices in space
 - ◇ Size of the slices or of the voxels
 - ◇ Slice offset — where is the dataset volume located in space?
 - ◇ For 3D+time datasets, also need slice timing information
 - ◇ to3d ‘knows’ how to get some of this auxiliary information from image file headers for some image file formats:
 - Mayo ANALYZE files (.hdr/.img pairs) contain voxel size information
 - Siemens .ima image arrays contain voxel size and orientation information
 - ◇ Other image formats either don’t have such information — or I just don’t know how to get it from the image files (yet)

- to3d runs in two modes:
 - ◇ Command line mode: you provide all auxiliary information on command line
 - Useful to get things done fast, and for AFNI experts
 - ◇ Graphical user interface (GUI) mode: you provide auxiliary information by filling out an on-screen form
 - Useful for explaining concepts, and for AFNI neophytes
- Sample #1: data from NIH GE 3 Tesla scanner
 - ◇ Files stored in archive AFNI_sample_01.tgz
 - unpack with command `gzip -dc AFNI_sample_01.tgz | tar xf -`
 - unpacks into directory AFNI_sample_01/
 - ◇ Anatomical (SPGR) data \implies 3D dataset (no time; 1 sub-brick)
 - 124 axial slices in subdirectory SPGR_anat
 - ◇ Functional (EPI) time series data \implies 3D+time dataset (160 sub-bricks)
 - 2880 images (18 coronal slices, 160 reps) in subdirectory EPI_run1
 - Visual stimulation task: rotating hemifield flashing checkerboard

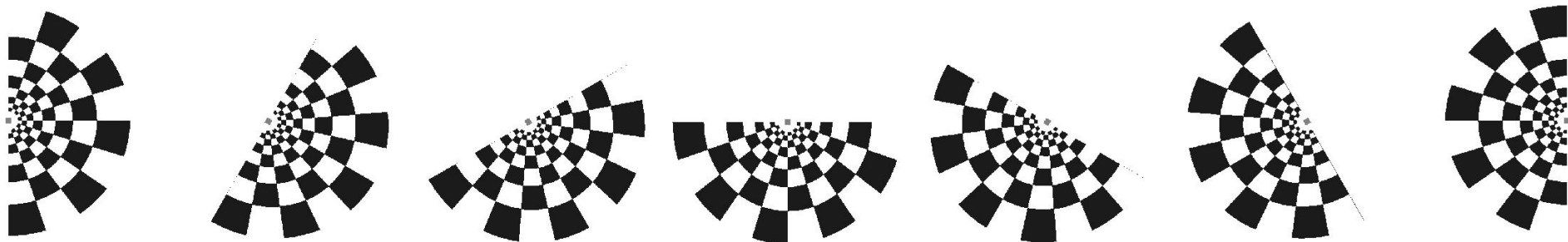

- Experiment log, taken at scanner:

NIMH-LBC-METHODS FMRI RESEARCH **RUN DATA**

ExpCode: AR25 Subject: [REDACTED] 30-18-08-8 Date: April 20 01 ~~PM~~ 8AM

Protocol#: F Start Time (24-hr): 8AM Investigator: Shruti Tech: ZSS

Scanner: 3T Coil: Wong / MAI / GE

Anat Scan 1 Type: (SPGR) / FSE / MPIR / IR / PSPGR TE(ms): M=full TR(ms): _____ Flip: 17° NEX: _____ FOV(mm): 240
 Matrix: 256 X 256 #Slices: 124 Plane: Ax / Cor / Sag Thickness(mm): 1.1 First: ISO.6 Last: 84.7 *Taken last*

Anat Scan 2 Type: SPGR / FSE / MPIR / TE(ms): _____ TR(ms): _____ Flip: _____ NEX: _____ FOV(mm): _____
 Matrix: X #Slices: _____ Plane: Ax / Cor / Sag Thickness(mm): _____ First: _____ Last: _____

EPI-Scan: (GE-EPI) / SE-EPI / GE RT EPI TE(ms): 25 TR(ms): 2 Flip: 90 NEX: _____ VTE: _____
 Plane: Ax / Cor / Sag FOV(mm): 240 Matrix: 64 Thickness(mm): 4 #Slices: 18 First: P96.8 Last: P28.8

Response Data: None / Accuracy / RT Other Data: EKG / HR / BP / RR / CO2 / GSR / EEG / EP / EMG / NONE _____

Timing: #Reps: _____ Skip#: _____ #Cycles: _____ #Conditions/Cycle: _____ #Reps/Condition/Cycle: _____
 Other Timing: 180 w/ slope 30 ON 15 OFF 270 total 30 pre 20 post ⇒ T.O.P

Run#	Time	Conditions	Stimulus File	Data File	Response File
		<u>180 w/ slope, CCW</u>	<u>003/I.019 → 043/I900</u>		
		<u>CCW</u>	<u>43/901 → 103/783</u>		
		<u>CW</u>	<u>103/784 → 163/666</u>		
		<u>CW</u>	<u>163/667 → 223/549</u>		<i>→ may have been ~ 1 sec starting scan ...</i>
		<u>CCW</u>	<u>223/550 → 283/432</u>		
		<u>CCW</u>	<u>283/433 → 343/315</u>		
		<u>CW</u>	<u>343/316 → 403/198</u>		
		<u>CW</u>	<u>403/199 → 463/081</u>		


(38s on command line, 8 sec. dropped)

- Using to3d to assemble the SPGR dataset [this is run on a Linux machine]:
 - ◇ `cd AFNI_sample_01/SGPR_anat` — change directory, to get at images
 - ◇ `ls` — to see what files are there (should see files I.001 ... I.124)
 - ◇ `to3d I.*` — run to3d, reading in all the image files — GUI pops up:




↪ N.B.: Warning about negative voxels appears for Linux/Intel computers, but not on Suns or SGIs

◇ To check images, click the View Images button in the to3d form:

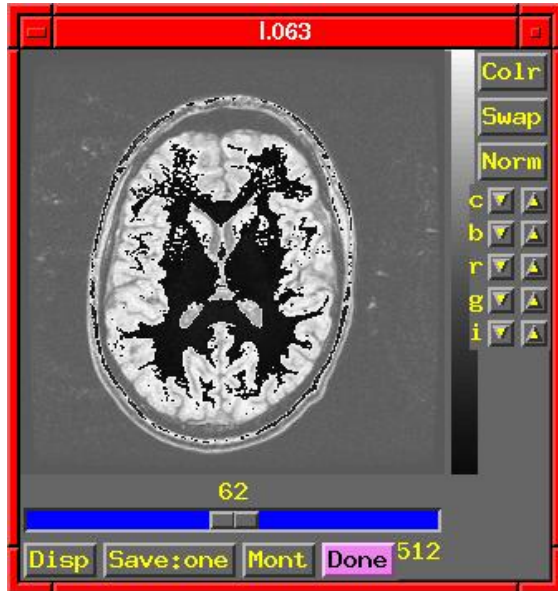

[On Linux/Intel computers]

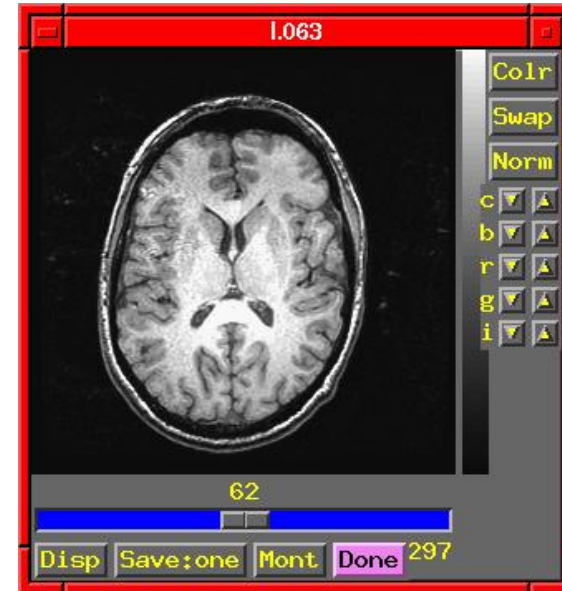

[On Sun or SGI computers]

◇ On Linux/Intel computers: the peculiar appearance of the image shows that something is wrong:

↳ MR images from scanners are stored as shorts: 2 bytes per number

↳ Like a 2 digit decimal number: "93" means $9 \times 10 + 3$

— By universal custom, we write the "9" first

— Could also write the same number as "39" (if we had a different custom)

↳ Customs for computers are not so universal

— Sun and SGI systems store 2 byte numbers in reverse order from Intel



— Result is that numbers are mangled (and some show up as negative)

— Solution: press to3d's Byte Swap[2] button, and images are fixed!


Same to3d control panel (without the negative voxel warning):



- ◇ Above the double line: must fill out 3 types of geometry information:
 - ↪ Left column: orientation of the dataset axes
 - ↪ Middle column: size of the dataset images or voxels
 - ↪ Right column: offset of the first slice

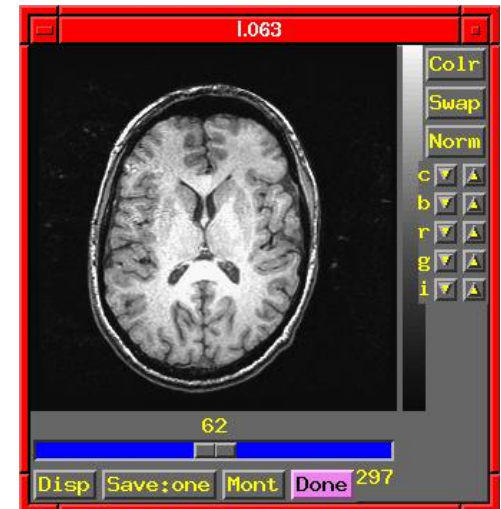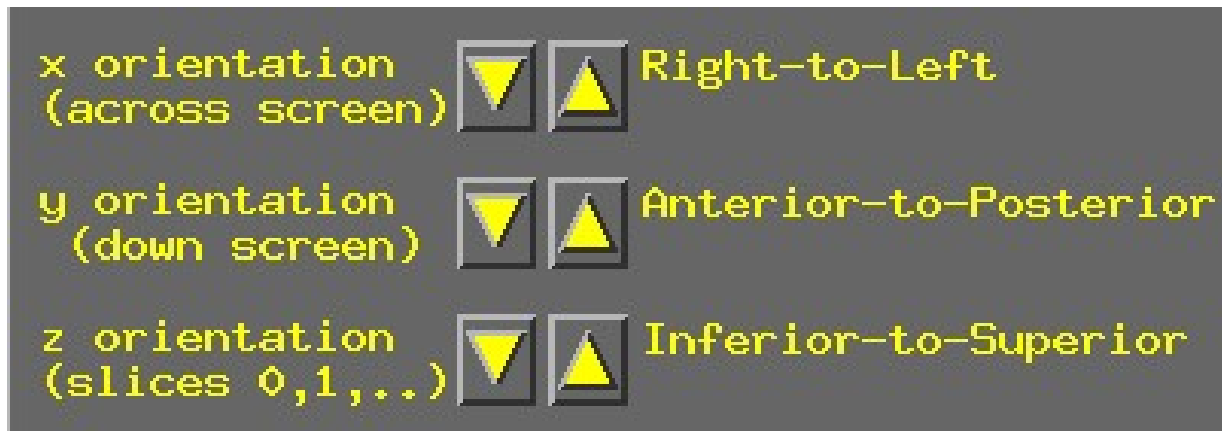

- ◇ Screen shot above shows correct orientation for this dataset
 - ↳ Use the image viewing window to judge how images are laid out
 - ↳ Click the arrows to scroll through the 6 possible options for each orientation to set the correct values
 - ↳ “x orientation” of dataset is across the screen (left to right)
 - item Must know subject’s right from left
 - ↳ “y orientation” of dataset is down the screen
 - ↳ “z orientation” of dataset is in increasing slice index order
 - determine this by using the slider at the bottom of image window

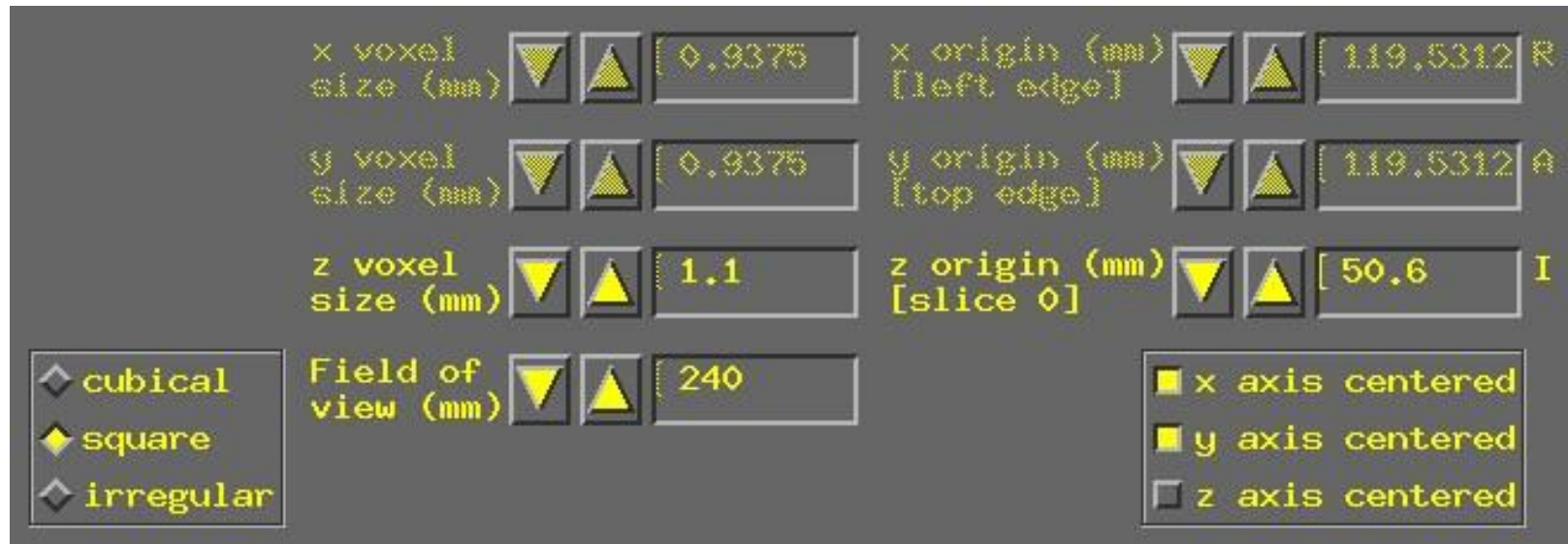

- ◇ To set dataset geometrical size/location, experiment log sheet is essential
- ◇ Screen shot above shows setting slice thickness to 1.1 mm
 - ↳ Default Field of view (FOV) of 240 mm is correct for these images
 - ↳ Default voxel geometry of “cubical” is incorrect
 - ↳ Must set geometry to “square” (x size = y size, z size varies)
 - ↳ Then set “z voxel size” to correct value (by typing in box)
- ◇ Screen shot shows setting location of center of first slice to 50.6 mm in Inferior (I) direction
 - ↳ Default is that slices are centered in the magnet
 - ↳ Probably not the case in the z direction
 - ↳ Click “z axis centered” off
 - ↳ Enter offset (here, 50.6 mm) into the “z origin” box


◇ Final required steps:

- ↪ Enter prefix for new dataset into “Prefix” text box at lower right of to3d control window
 - Choosing a good prefix is important for keeping datasets organized
- ↪ Press “Save Dataset” button
- ↪ Press “quit” (twice) to exit to3d
- ↪ The new dataset files should show up when you use command ls
- ↪ You might want to move them to some other directory
 - mv *+orig.* ../afni to move datasets to directory named afni, one level above
 - this directory was created when you unpacked AFNI_sample_01.tgz; it contains pre-made AFNI datasets from EPI and SPGR images

◇ Other image sources (besides GE reconstruction):

- ↪ Mayo ANALYZE files (.hdr/.img pairs) contain voxel size information (orientation information not always reliable)
- ↪ Siemens .ima image arrays contain voxel size and orientation information
- ↪ Some day I may learn how to read such information from GE image headers!
- ↪ Can also assemble datasets of bytes and floats
- ↪ If all else fails when reading an image, see AFNI FAQ #66

- Using to3d to assemble the EPI dataset
 - ◇ cd ../EPI_run1 — change directory, to get at images
 - ◇ ls — to see what files are there (should see files I.00001 ... I.02880)
 - ◇ We do not just do to3d I.* to create a 3D+time dataset
 - ◇ For hysterical historical reasons, the time-axis information must be given on the to3d command line
 - Cannot be modified from GUI
 - ◇ to3d -time:zt 18 160 2000 alt+z I.*
 - ◇ -time:zt means slices will be presented in order of space (z) then time (t)
 - This is the usual way slices are ordered, but -time:tz is needed at some sites
 - If in doubt, do to3d I.*, use viewer to look at slices and see their order
 - ◇ 18 160 means that there will be 18 slices in z, 160 in t (2880 total)
 - ◇ 2000 means that the TR for volume acquisition was 2000 ms
 - Could also do 2s to specify that TR is in seconds
 - ◇ alt+z means that the slices are gathered in alternating order in the +z direction
 - Most EPI acquisitions are really 2D multislice, spread out through time
 - AFNI header contains information about slice timing offsets
 - Other possible modes: zero (for 3D), @filename (to specify each slice)

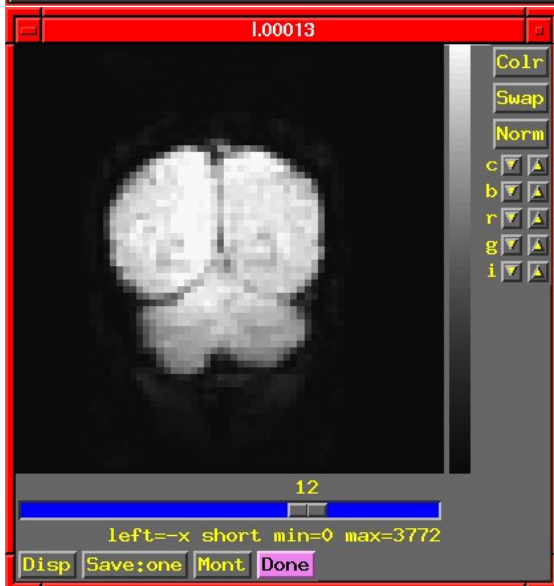

- ◇ Fields in GUI screenshot above are already filled out — Coronal slices; will work with SPGR axials in AFNI
- ◇ Note slice thickness and slice offset (“z origin”) — Taken from experiment log
- ◇ Time information is displayed, but not editable
- ◇ Have set “Type of anatomy” to “Echo Planar” — Just acts as a reminder to user